

# Bisulfite Converted Universal Methylated Human DNA Standard

Cat. Nos. D5015

Storage: -20 °C



## Product Information

### Highlights:


- Purified male human DNA that has been methylated at all cytosine positions within the CG dinucleotide context by M.SssI methyltransferase and treated with sodium bisulfite for direct downstream methylation applications.
- Control primer pair is designed to amplify a fragment of the human MLH1 mismatch repair gene following bisulfite conversion.

### Product Contents:

	Cat. # D5015	Storage Temp.
Bisulfite Converted Universal Methylated Human DNA Standard	1 µg/50 µl	-20 °C
hMLH1 Primers	20 µl	-20 °C

### Description:

The **Bisulfite Converted Universal Methylated Human DNA Standard** contains bisulfite-treated DNA derived from the **Universal Methylated Human DNA Standard**. The supplied DNA has been methylated at all cytosine positions within the CG dinucleotide context using M.SssI methyltransferase<sup>1</sup> (EC 2.1.1.37; Figure 1). Following methylation, the DNA was treated with sodium bisulfite according to the protocol described in our **EZ DNA Methylation-Direct™** kit. In this kit the methylated cytosines are unmodified by sodium bisulfite, whereas all non-methylated cytosines are converted to uracil, which are detected as thymine in subsequent PCR.

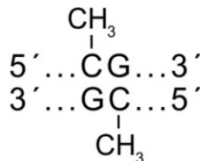

**Figure 1.** M.SssI methyltransferase methylates all cytosine residues in the double-stranded CpG context.

The included primer set is designed to amplify a fragment of the human MLH1 mismatch repair gene (GenBank Accession# U83845, nucleotides 804 through 986) and can be used directly with the Bisulfite Converted Universal Methylated Human DNA Standard.

### References:

- Nur *et al.* J. Bacteriol. 164: 19-24 (1985).

### Protocol:

*Note: We recommend using ZymoTaq™ DNA Polymerase or other hot-start DNA polymerases for amplification of bisulfite-treated DNA.*

#### 1. PCR Setup:

The following setup is designed for a 25 µl total reaction volume:

Component	Volume	Final Conc.
hMLH1 primer I*	Variable	0.2 to 0.8 µM
hMLH1 primer II*	Variable	0.2 to 0.8 µM
Bisulfite-converted DNA	1 µl	0.8 ng/µl
10 mM dNTP mix	0.5 µl	0.2 mM each dNTP
Standard PCR buffer	Variable	1x
MgCl <sub>2</sub> or MgSO <sub>4</sub>	Variable	1-4 mM, if needed
ZymoTaq™ DNA Polymerase (or other Hot-start DNA polymerase)	Variable	1 to 2 units
Add water to 25 µl		

\* Alternatively, you may substitute primers of your choice.

#### 2. Recommended Thermocycler Conditions:

- 95 °C, 10 minutes
- 95 °C, 30 seconds
- 59 °C, 30 to 60 seconds
- 72 °C, 60 seconds
- Repeat steps B through D an additional 29 to 39 times depending on the polymerase used.
- 72 °C, 7 minutes
- 4 °C

The PCR amplicon can now be used directly for sequencing analysis or cloning.

### Product Specifications:

- Universal Methylated DNA Standard, 1 µg/50 µl.  
Source: DNA isolated from male human enzymatically methylated by M.SssI Methyltransferase (EC 2.1.1.37) and bisulfite converted.  
Concentration: 20 ng/µl in TE buffer (10 mM Tris-HCl, 0.1 mM EDTA, pH 8.0).  
Storage: -20 °C
- Control Primers.  
Concentration: 20 µM each primer in TE buffer (10 mM Tris-HCl, 1 mM EDTA, pH 8.0)  
Volume: 20 µl of mixed primers  
Storage: -20 °C  
Sequence:

hMLH1 Primer I:

5' - GGAGTGAAGGAGGTTACGGGTAAGT - 3'

hMLH1 Primer II:

5' - AAAAACGATAAAACCCTATACCTAATCTATC - 3'

Continued on next page...

## Appendix:

The expected PCR amplicon for the Bisulfite Converted Universal Methylated Human DNA Standard is 182 bp, corresponding to nucleotide positions 804 to 986 of human MLH1 DNA including the regions (italicized) that hybridize to the primers (GenBank Accession #: U83845).

Methylated cytosines (underlined> in the CpG dinucleotide context (bold capital letters) remain unconverted following bisulfite treatment, whereas non-methylated cytosines, or cytosines not in the CpG context, are converted to uracil and detected as thymine after PCR.

```

801  ---ggagtga aggaggttaC GggttaagtCG ttttgaCGta
841  gaCGttttat tagggtCGCG CGttCGtCGt tCGttatata
881  tCGttCGtag tattCGtggt tagtttCGta gtggCGtttg
921  aCGtCGCGtt CGCGgtagt taCGatgagg CGCGataga
961  ttaggtatag ggttttatCG ttttt-----

```

## Trademarks and Disclaimers:

™ Trademarks of Zymo Research Corporation.

This product is for research use only and should only be used by trained professionals. Wear protective gloves and eye protection. Follow the safety guidelines and rules enacted by your research institution or facility.

The Polymerase Chain Reaction (PCR) process is covered by U.S. Patent: #4,683,195; 4,683,202 assigned to Hoffmann-La Roche. Patents pending in other countries. No license under these patents to use the PCR process is conveyed expressly or by implication to the purchaser by the purchase of Zymo Research's products. Further information on purchasing licenses to practice the PCR process can be obtained from the director of Licensing at Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404 or at Roche Molecular Systems, Inc., 1145 Atlantic Avenue, Alameda, California 94501.

Version 1.0.5

## Also Available:

Product Name	Size	Catalog number
EZ DNA Methylation™ Kit	50	D5001
	200	D5002
	2 x 96	D5003
	2 x 96	D5004
EZ DNA Methylation-Gold™ Kit	50	D5005
	200	D5006
	2 x 96	D5007
	2 x 96	D5008
EZ DNA Methylation-Direct™ Kit	50	D5020
	200	D5021
	2 x 96	D5022
EZ DNA Methylation-Startup™ Kit	2 x 96	D5023
	1 Kit	D5024
EZ Bisulfite DNA Clean-up Kit™	50	D5025
	200	D5026
	2 x 96	D5027
	2 x 96	D5028
Universal Methylated DNA Standard	1 set	D5010
Universal Methylated Human DNA Standard	1 set	D5011
Universal Methylated Mouse DNA Standard	1 set	D5012
Human HCT116 DKO Methylation Standards	1 set	D5014
Human HCT116 DKO Non-methylated DNA Standard	5 µg	D5014-1
Human HCT116 DKO Methylated DNA Standard	5 µg	D5014-2
<i>E. coli</i> Non-methylated Genomic DNA	5 µg	D5016
ChIP DNA Clean & Concentrator™	50	D5201
	50	D5205
Methylated-DNA IP Kit	10	D5101
Anti-5-Methylcytosine Monoclonal Antibody (clone 10G4)	50 µg 200 µg	A3001-50 A3001-200
ZymoTaq™ DNA Polymerase	50	E2001
	200	E2002
ZymoTaq™ PreMix (2X concentrated)	50	E2003
	200	E2004
CpG Methylase (M.SssI)	200 units	E2010
	400 units	E2011

## ZYMO RESEARCH CORP.

Phone: (949) 679-1190 ▪ Toll Free: (888) 882-9682 ▪ Fax: (949) 266-9452 ▪ info@zymoresearch.com ▪ www.zymoresearch.com