



ZYMO RESEARCH

The Beauty of Science is to Make Things Simple

INSTRUCTION MANUAL

EZ Bisulfite DNA Clean-up Kit™

Catalog Nos. **D5025 & D5026**

Highlights

- Consistent, high yield recovery of bisulfite-treated DNA from any “homebrew” or commercial reaction mixture containing bisulfite.
- Simple spin column procedure with small elution volumes ($\geq 10 \mu\text{l}$) for concentrated DNA.
- Recovered DNA is ideal for downstream methylation analysis including PCR, endonuclease digestion, sequencing, arrays, etc.

Contents

Product Contents	1
Product Specifications.....	1
Introduction to DNA Methylation.....	2
Product Description.....	3
Reagent Preparation	4
Protocol.....	4
Appendix	5
Frequently Asked Questions	6
Ordering Information	7
List of Related Products	8

For Research Use Only

Ver. 1.0.1

ZYMO RESEARCH CORP.

Phone: (949) 679-1190 ▪ Toll Free: (888) 882-9682 ▪ Fax: (949) 266-9452 ▪ info@zymoresearch.com ▪ www.zymoresearch.com

Product Contents:

Note: Satisfaction of all Zymo Research products is guaranteed. If you should be dissatisfied with this product please call 1-888-882-9682.

EZ Bisulfite DNA Clean-up Kit™ (Format)	D5025 (50 preps.)	D5026 (200 preps.)	Storage Temperature
M-Binding Buffer	20 ml	80 ml	Room Temp.
M-Wash Buffer¹	6 ml	24 ml	Room Temp.
M-Desulphonation Buffer	10 ml	40 ml	Room Temp.
M-Elution Buffer	1 ml	4 ml	Room Temp.
Zymo-Spin™ IC Columns	50 columns	200 columns	Room Temp.
Collection Tubes	50 tubes	200 tubes	Room Temp.
Instruction Manual	1	1	–

Note - Integrity of kit components is guaranteed for one year from date of purchase. Reagents are routinely tested on a lot-to-lot basis to ensure they provide maximal performance and reliability.

¹ Add 24 ml of 100% ethanol to the 6 ml **M-Wash Buffer** concentrate (D5025) or 96 ml of 100% ethanol to the 24 ml **M-Wash Buffer** concentrate (D5026) before use.

Specifications:

- **DNA Input:** Up to 5 µg DNA (per column) from “homebrew” or commercial reaction mixtures containing bisulfite.
- **DNA Recovery:** Typical yields are between 80-90%.
- **Elution:** DNA can be eluted with ≥ 10 µl M-Elution Buffer or water.
- **Application:** Purified DNA can be used directly for PCR, endonuclease digestion, arrays, sequencing, etc.

The Polymerase Chain Reaction (PCR) process is covered by U.S. Pat. Nos. 4,683,195 and 4,683,202 assigned to Hoffmann-La Roche. Patents pending in other countries. No license under these patents to use the PCR process is conveyed expressly or by implication to the purchaser by the purchase of Zymo Research's EZ DNA Methylation kits. Further information on purchasing licenses to practice the PCR process can be obtained from the director of Licensing at Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404 or at Roche Molecular Systems, Inc., 1145 Atlantic Avenue, Alameda, California 94501.

Use of Methylation Specific PCR (MSP) is protected by US Patents 5,786,146 & 6,017,704 & 6,200,756 & 6,265,171 and International Patent WO 97/46705. No license under these patents to use the MSP process is conveyed expressly or by implication to the purchaser by the purchase of this product.


Note - ™ Trademarks of Zymo Research Corporation. This product is for research use only and should only be used by trained professionals. Some reagents included with this kit are irritants. Wear protective gloves and eye protection. Follow the safety guidelines and rules enacted by your research institution or facility.

ZYMO RESEARCH CORP.

Introduction to DNA Methylation:

DNA methylation is a naturally occurring event in both prokaryotic and eukaryotic organisms. In prokaryotes DNA methylation provides a way to protect host DNA from digestion by restriction endonucleases that are designed to eliminate foreign DNA, and in higher eukaryotes DNA methylation functions in the regulation/control of gene expression (1). It has been demonstrated that aberrant DNA methylation is a widespread phenomenon in cancer and may be among the earliest changes to occur during oncogenesis (2). DNA methylation has also been shown to play a central role in gene imprinting, embryonic development, X-chromosome gene silencing, and cell cycle regulation. In many plants and animals, DNA methylation consists of the addition of a methyl group to the fifth carbon position of the cytosine pyrimidine ring via a methyltransferase enzyme (3). The majority of DNA methylation in mammals occurs in 5'-CpG-3' dinucleotides, but other methylation patterns do exist. In fact, about 80 percent of all 5'-CpG-3' dinucleotides in mammalian genomes are found to be methylated, whereas the majority of the twenty percent that remain unmethylated are within promoters or in the first exons of genes.

The ability to detect and quantify DNA methylation efficiently and accurately has become essential for the study of cancer, gene expression, genetic diseases, as well as many other important aspects of biology. To date, a number of methods have been developed to detect/quantify DNA methylation including: high-performance capillary electrophoresis (4) and methylation-sensitive arbitrarily primed PCR (5). However, the most common technique used today remains the bisulfite conversion method (6). This technique involves treating methylated DNA with bisulfite, which converts unmethylated cytosines into uracil. Methylated cytosines remain unchanged during the treatment. Once converted, the methylation profile of the DNA can be determined by PCR amplification followed by DNA sequencing (see below) or other analytical procedures.

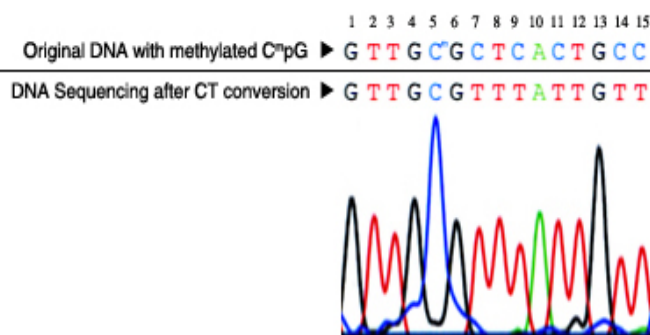

DNA sequencing results following bisulfite treatment. DNA with methylated C^mpG at nucleotide position #5 was processed using the **EZ DNA Methylation™ Kit**. The recovered DNA was amplified by PCR and then sequenced directly. The methylated cytosine at position #5 remained intact while the unmethylated cytosines at positions #7, 9, 11, 14 and 15 were completely converted into uracil following bisulfite treatment and detected as thymine following PCR.

References:

1. Costello JF, Plass CJ. *Med. Genet.* 2001; 38(5): 285-303.
2. Stirzaker C. *Cancer Res.* 1997; 57(11): 2229-2237.
3. Adams RL. *Bioessays.* 1995; 17(2): 139-145.
4. Fraga MF, *et al.* *Electrophoresis.* 2000; 21(14): 2990-2994.
5. Gonzalgo ML. *Cancer Res.* 1997; 57(4): 594-599.
6. Frommer M. *Proc. Natl. Acad. Sci. USA.* 1992; 89(5): 1827-1831.

Selected EZ DNA

Methylation™ Kit Citations:

1. Ehrich M, *et al.* Nuc. Acids Res. 2007; 35 (5): e29
2. Kaneda M, *et al.* Nature. 2004; 429: 900-903
3. Zhang F, *et al.* Proc. Natl. Acad. Sci. USA. 2007; 104 (11): 4395-4400.
4. Oda M, *et al.* Genes & Dev. 2006; 20: 3382-3394.
5. England RPM, *et al.* Nature Meth. 2005; 2: 1-2.

Product Description:

The **EZ Bisulfite DNA Clean-up Kit™** has been specifically designed for the purification of bisulfite-treated DNA from any “homebrew” or commercial reaction mixture containing bisulfite.

The product features innovative *Fast-Spin*, in-column desulphonation and wash technologies that eliminate DNA loss, buffer carryover and the need for ethanol/isopropanol precipitations. The procedure is easy and DNA clean-up can be completed in just minutes. Bisulfite-treated DNA purified with the **EZ Bisulfite DNA Clean-up Kit™** is ideal for PCR amplification for downstream DNA methylation analysis including endonuclease digestion, sequencing, microarrays, etc.



DNA is perfect for...

- PCR
- Endonuclease Digestion
- Arrays
- Sequencing

Buffer Preparation:

Add 24 ml of 100% ethanol to the 6 ml **M-Wash Buffer** concentrate (D5025) or 96 ml of 100% ethanol to the 24 ml **M-Wash Buffer** concentrate (D5026) before use.

Protocol:

1. Add 4 volumes of **M-Binding Buffer** to each volume of a bisulfite-containing reaction mixture (4:1) and mix.

For Example: Add 400 μ l M-Binding Buffer to 100 μ l from a conversion reaction.

2. Load the mixture into a **Zymo-Spin™ IC Column** in a **Collection Tube**.
3. Centrifuge at 10,000 x *g* for 30 seconds. Discard the flow-through.
4. Add 100 μ l of **M-Wash Buffer** to the column. Centrifuge at 10,000 x *g* for 30 seconds.
5. Add 200 μ l of **M-Desulphonation Buffer** to the column and let stand at room temperature (20°C - 30°C) for 15 - 20 minutes. After the incubation, centrifuge at 10,000 x *g* for 30 seconds.
6. Add 200 μ l of **M-Wash Buffer** to the column. Centrifuge at 10,000 x *g* for 30 seconds. Add another 200 μ l of **M-Wash Buffer** and centrifuge for an additional 30 seconds.
7. Place the column into a 1.5 ml microcentrifuge tube. Add 10 μ l of **M-Elution Buffer** directly to the column matrix. Centrifuge for 30 seconds at 10,000 x *g* to elute the DNA.

The capacity of the collection tube with the column inserted is 800 μ l. Empty the collection tube whenever necessary to prevent contamination of the column contents by the flow-through.

Alternatively, water or TE (pH \geq 6.0) can be used for elution if required for your experiments.

The DNA is ready for immediate analysis or can be stored at or below -20°C for later use. For long term storage, store at or below -70°C. The elution volume can be > 10 μ l depending on the requirements of your experiments, but small elution volumes will yield more concentrated DNA.

Important: PCR of bisulfite-treated DNA is from single stranded templates, since most of the DNA is single stranded and non-complementary following bisulfite treatment.

Zymo Taq™ is a "hot start" DNA polymerase specifically designed for the amplification of bisulfite treated DNA. (see page 8 for details)

Appendix: PCR of Bisulfite-Treated DNA


1. **PCR Primer Design.** Generally, primers of 24 to 32 bases are required for amplification of bisulfite converted DNA. For most eukaryotes, all non-methylated cytosine residues will be converted into uracil during the bisulfite treatment. These Cs should be treated as Ts for primer design purposes. For example, for the sequence 5'-AACCTTACAGGCAC-3', the corresponding primer should be 5'-AA**TTTTATAGGTAT**-3'.

If the primer contains CpG dinucleotides with uncertain methylation status, then mixed bases with C and T can be used. Usually, there should be no more than three mixed positions per primer and they should be located toward the 5' end of the primer. It is not recommended to have mixed bases located at the 3' end of the primer.

2. **PCR Conditions.** Usually, 35 to 40 cycles are required for successful PCR amplification of bisulfite converted DNA. Optimal amplicon size should be between 150 - 300 bp; however larger amplicons (up to 1 kb) can be generated with optimization of the bisulfite reaction and PCR conditions. We have found that annealing temperatures between 55 - 60°C typically work well. As most non-methylated cytosine residues are converted into uracil, the bisulfite-treated DNA usually is AT-rich and has low GC composition. Thus, it may be necessary to reduce the annealing temperature accordingly.

Non-specific PCR amplification is relatively common with bisulfite treated DNA due to its AT-rich nature. PCR using "hot start" polymerases is strongly recommended for the amplification of bisulfite-treated DNA.

3. **Quantifying Bisulfite-Treated DNA.** Following bisulfite treatment of genomic DNA, non-methylated cytosine residues are converted into uracil. The recovered DNA is typically A, U, and T-rich. The original base-pairing no longer exists. Instead, it is single stranded with limited non-specific base-pairing at room temperature. The absorption coefficient at 260 nm resembles that of RNA. Use a value of 40 µg/ml for $A_{260} = 1.0$ when determining the concentration of the recovered bisulfite-treated DNA.
4. **"Visualizing" Bisulfite-Treated DNA.** Bisulfite-treated DNA can be visualized in agarose/EtBr gels following electrophoresis using a standard UV-light source. However, cooling the gel on ice for 10-15 minutes prior to visualization will greatly enhance the resolution and apparent banding of the DNA (see below).

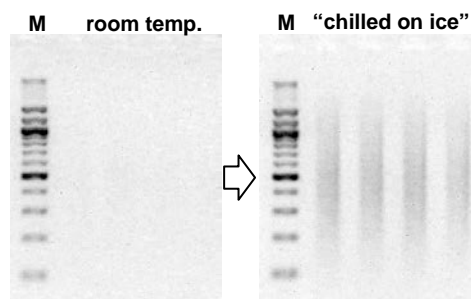

"Visualizing" bisulfite-treated DNA in agarose/EtBr gels is best done following chilling the gels on ice (see 4. on page 5). In the figures above, bisulfite-treated salmon sperm DNA was desulphonated then purified using the **EZ Bisulfite DNA Clean-up Kit™**. The DNA, mostly single stranded, was then separated in a 0.8 % (w/v) agarose/TAE/EtBr gel and visualized with a UV-light source immediately following electrophoresis (room temp) and after chilling the gel on ice for 15 minutes. M is a 100 bp DNA ladder (Zymo Research).

Frequently Asked Questions:

Q: At what temperature and for how long can converted DNA be stored?


A: *The sample should be stored at $\leq -20^{\circ}\text{C}$ whenever possible. The quality of the DNA should remain relatively unchanged for up to 3 months.*

Q: Which *Taq* polymerase(s) do you recommend for PCR amplification of converted DNA?

A: *We recommend a “hot start” DNA polymerase (e.g., ZymoTaq™, page 8).*

Q: Why are there two different catalog numbers for the EZ-96 Bisufite DNA Clean-up Kit™?

A: *The two different catalog numbers are used to differentiate between the binding plates that are included in the kit. Deep and shallow-well binding plates are available to accommodate most rotors and microplate carriers. Below is a comparison of the two binding plates.*

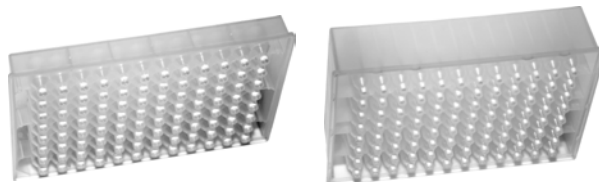

Binding Plate	Silicon-A™ Plate	Zymo-Spin™ I-96 Plate
Style	Shallow-Well	Deep-Well
Height of Binding Plate	19 mm (0.75 inches)	35 mm (1.38 inches)
Binding Plate/Collection Plate Assembly	43 mm (1.69 inches)	60 mm (2.36 inches)
Binding Cap./Minimum Elution Volume	5 µg/30 µl	5 µg/15 µl
Catalog Numbers	D5027	D5028

Ordering Information:

Product Description	Catalog No.	Kit Size
EZ Bisufite DNA Clean-up Kit™	D5025	50 preps.
	D5026	200 preps.
EZ-96 Bisufite DNA Clean-up Kit™ (Shallow-Well)	D5027	2 x 96 preps.
EZ-96 Bisufite DNA Clean-up Kit™ (Deep-Well)	D5028	2 x 96 preps.

For Individual Sale	Catalog No.	Amount(s)
M-Binding Buffer	D5001-3	20 ml
	D5002-3	80 ml
M-Wash Buffer	D5001-4	6 ml
	D5002-4	24 ml
	D5007-4	36 ml
M-Desulphonation Buffer	D5001-5	10 ml
	D5002-5	40 ml
M-Elution Buffer	D5001-6	1 ml
	D5002-6	4 ml
Zymo-Spin™ IC Columns (capped)	C1004-50	50 columns
	C1004-250	250 columns
Collection Tubes	C1001-50	50 tubes
	C1001-500	500 tubes
	C1001-1000	1,000 tubes

ZYMO RESEARCH CORP.

Epigenetics Products From Zymo Research

Product	Description	Kit Size	Cat No. (Format)
Bisulfite Kits for DNA Methylation Detection			
EZ DNA Methylation™ Kit	For the conversion of unmethylated cytosines in DNA to uracil via the <u>chemical-denaturation</u> of DNA and a specially designed CT Conversion Reagent. <i>Fast-Spin</i> technology ensures ultra-pure, converted DNA for subsequent DNA methylation analysis.	50 Rxns. 200 Rxns. 2x96 Rxns. 2x96 Rxns.	D5001 (spin column) D5002 (spin column) D5003 (shallow-well plate) D5004 (deep-well plate)
EZ DNA Methylation-Gold™ Kit	For the fast (3 hr.) conversion of unmethylated cytosines in DNA to uracil via <u>heat/chemical-denaturation</u> of DNA and a specially designed CT Conversion Reagent. <i>Fast-Spin</i> technology ensures ultra-pure, converted DNA for subsequent DNA methylation analysis.	50 Rxns. 200 Rxns. 2x96 Rxns. 2x96 Rxns.	D5005 (spin column) D5006 (spin column) D5007 (shallow-well plate) D5008 (deep-well plate)
EZ DNA Methylation-Direct™ Kit	Features simple and reliable DNA bisulfite conversion directly from blood, tissue (FFPE/LCM), and cells without the prerequisite for DNA purification in as little as 4-6 hrs. The increased sensitivity of this kit makes it possible to amplify bisulfite converted DNA from as few as 10 cells or 50 pg DNA.	50 Rxns. 200 Rxns. 2x96 Rxns. 2x96 Rxns.	D5020 (spin column) D5021 (spin column) D5022 (shallow-well plate) D5023 (deep-well plate)
EZ DNA Methylation-Startup™ Kit	Designed for the first time user requiring a consolidated product to perform DNA methylation analysis. Includes technologies for sample processing, bisulfite treatment of DNA, and PCR amplification of "converted" DNA for methylation analysis.	1 Kit	D5024
EZ Bisulfite DNA Clean-up Kit™	Desulphonation and purification of DNA from any "homebrew" or commercially derived reaction mixture containing bisulfite.	50 Preps. 200 Preps. 2x96 Preps. 2x96 Preps.	D5025 (spin column) D5026 (spin column) D5027 (shallow-well plate) D5028 (deep-well plate)
Methylated DNA Standards			
Universal Methylated DNA Standard	pUC19 plasmid DNA having all CpG sites methylated. To be used for the evaluation of bisulfite-mediated conversion of DNA. Supplied with a control primer set.	1 set	D5010
Universal Methylated Human DNA Standard	Human (male) genomic DNA having all CpG sites methylated. To be used for the evaluation of bisulfite-mediated conversion of DNA. Supplied with a control primer set.	1 set	D5011
Universal Methylated Mouse DNA Standard	Mouse (male) DNA having all CpG sites methylated. To be used for the evaluation of bisulfite-mediated conversion of DNA. Supplied with a control primer set.	1 set	D5012
Other...			
ChIP DNA Clean & Concentrator™	Clean and concentrate DNA from any reaction or "crude" preparation in 2 min. A 6 µl minimum elution volume allows for highly concentrated DNA. Designed for samples containing up to 5 µg of DNA.	50 Preps. 50 Preps.	D5201 (uncapped column) D5205 (capped column)
ZymoTaq™ DNA Polymerase	ZymoTaq™ "hot start" DNA Polymerase is specifically designed for the amplification of "difficult" DNA templates including: bisulfite-treated DNA for methylation detection. The product generates specific amplicons with little or no by-product formation. Available either as a single buffer premix or as a polymerase system with components provided separately.	50 Rxns. 200 Rxns. 50 Rxns. 200 Rxns.	E2001 (system) E2002 (system) E2003 (premix) E2004 (premix)
Anti-5-Methylcytosine Monoclonal Antibody (clone 10G4)	Mouse monoclonal antibody developed to facilitate the differentiation between methylated and non-methylated cytosines in DNA. Can be used in immunoprecipitation-based procedures including Methylated DNA Immunoprecipitation (MeIP).	50 µg/50 µl 200 µg/200 µl	A3001-50 A3001-200
Methylated-DNA IP Kit	IP with a highly specific anti-5-methylcytosine monoclonal antibody. Designed for the enrichment of 5-methylcytosine-containing DNA from any pool of fragmented genomic DNA for use in genome-wide methylation analysis.	10 Rxns.	D5101

ZYMO RESEARCH CORP.